



Evolve Scope Drying Cabinet with Hinged Door

Evolve modular cabinets and casework create a flexible storage system that organizes a hospital's most supply-intensive departments. A clean design and tight construction offer a refined aesthetic and support infection prevention efforts.

The 19.25"-deep scope drying cabinet has a tempered glass door and comes in 2 widths. The cabinet's drying mechanism continuously filters air into all internal channels to more quickly and thoroughly dry scopes; the HEPA filter offers a 99.97% retention rate of particles less than 0.3µm. It has an automated drying cycle, self-adjusting high-pressure regulation, automatic high-pressure shut-off, and fulfills forced air recommendations. The cabinet material is a nonporous AireCore with Glacier White finish or with a specified laminate surface. Cabinet has a 4"-high toe kick.

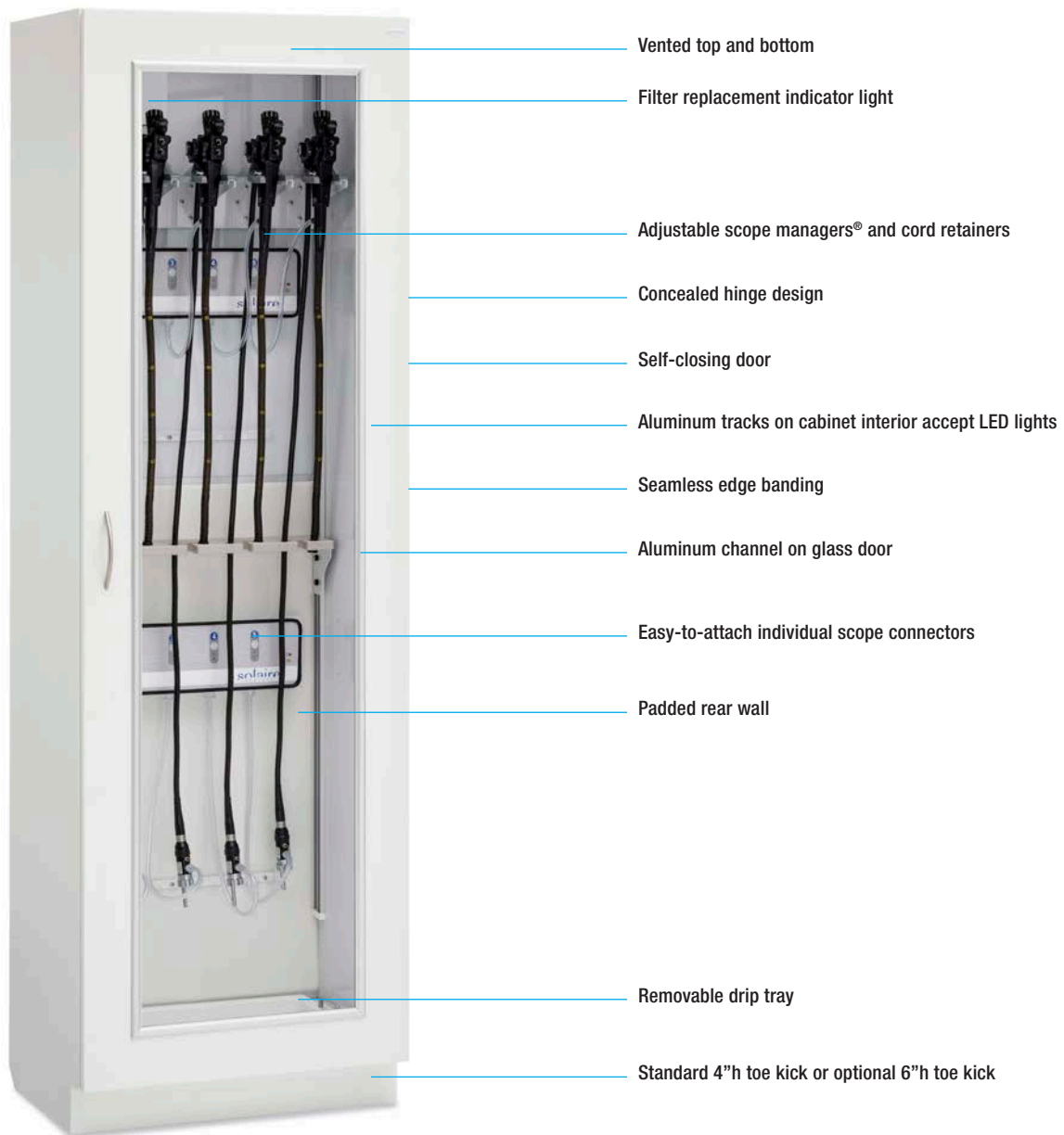

Evolve Scope Drying Cabinet with Hinged Door

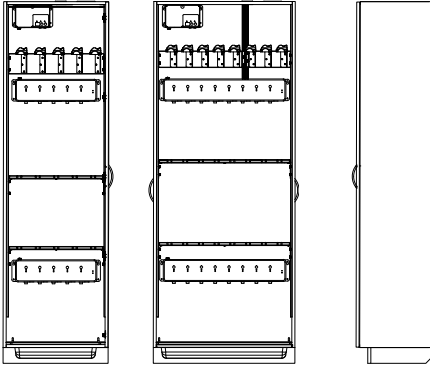

Dimensions (Dimensions are nominal)

Exterior (interior); scope hanging length:

26.75"w x 19.25"d x 92"h (25.25"w x 17"d x 86.5"h); 83" maximum (shown)

36"w x 19.25"d x 92"h (34.5"w x 17"d x 86.5"h); 83" maximum (shown)

92"h cabinet with 6"h toe kick is 94"h

Interior Configuration:

92"h x 26.75"w cabinet holds 5 or 10 scopes

92"h x 36"w cabinet holds 8 or 16 scopes

Options

- Door can be hinged on right or left side of 26.75"w cabinet
- Optional Fan Ventilation Unit available
- Keyed or keyless lock options available
- Optional sloped top available
- Optional 6" toe kick



Product Numbers

26.75"w x 19.25"d x 92"h

AireCore: SEA1927SGLSF, SEA1927SGRSF (shown)

With double scope manager: SEA1927SGLDF, SEA1927SGRDF

36"w x 19.25"d x 92"h

AireCore: SEA1936SGSF

With double scope manager: SEA1936SGDF

Materials and Finishes

Cabinet shell: **Aluminum composite with Glacier White finish or specified Wilsonart, Formica, or Pionite standard laminate**

Edge banding: **White ABS**

Glass door: **Clear anodized aluminum and tempered glass**

Glass trim: **Clear anodized aluminum**

Hinges: **Metal concealed European-style**

Door pulls: **Dull chrome-plated steel**



Storage solutions for healthcare

Solairé Medical | P.O. Box 2808 | Grand Rapids, Michigan 49501
888-435-2256 phone | 616-743-5944 fax | solairomedical.com